

Enzo Pendant Technical Guide



Introduction

Our most popular personal alarm pendant, Enzo, is a lightweight device that can be worn on the wrist or around the neck. It has a long battery life and is waterproof and shockproof.

Enzo has two-way data communication which means that the user receives a confirmation that the alarm has been received at our contact centre. When the button is pressed, the LED indicator lights up. When the alarm base unit receives the signal from the button, this is confirmed by the indicator light changing color.

Our unique Easy Press Concept (EPC) accessory makes it easy for users with restricted movement to call for assistance.

Enzo is compatible with all of Careium's current alarm base units: CareIP®Mobile, CareMobile, Eliza & Eliza S.

Description of functions

The Enzo personal alarm pendant can be worn around the wrist or around the neck. The Easy press accessory works in both positions. You can easily switch between the different accessories.

- LED indicates transmission and battery status.
- Range covers a normal house or flat.
- Confirmation of successful alarm transmission.
- Lithium battery for long life, up to 5 years.
- Automatic monitoring of battery and function. Careium is alerted in the event of a fault.
- Battery replaceable by authorised personnel.
- Adjustable neck strap.



Technical data

ENZO PERSONAL ALARM PENDANT	
Frequency	869.2 MHz
Water resistance	According to IP67 standard
Battery	1 Lithium, type CR 2032
Automatic alarm test	Yes, selectable
Dimensions	39 x 32 x 11 mm
Article number	300 208

Careium Enzo is also available in a variant for fixed wall mounting. For example, in the shower or on the night table.

Article number 300212.

Approval



-marked. Complies with relevant standards.

EU declaration of conformity

Hereby, Careium declares that the radio equipment type is in compliance with Directives: 2014/53/EU and 2011/65/EU. The full text of the EU declaration of conformity is available at the following internet address: www.careium.com/dofc

Product information

This product is a part of a social alarm system intended to initiate a call in case of emergency or if assistance is needed.

- Frequency band: 869.200-869.250 MHz
- Radio-frequency power: max. 10 mW (10 dBm)

Safety information


- The unit is waterproof for use in shower or bath (IP67).
- Battery replacement may be performed only by authorised persons.

Environmental information


This product complies with the requirements of the EU directive 2006/66/EC (Batteries) and 2012/19/ EU (WEEE) and should not be disposed of with other household waste.

Always recycle your used electronic products, batteries, and packaging materials to dedicated collection points.



1

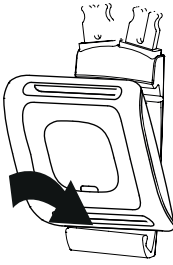


2

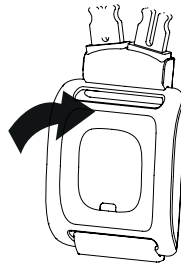


3

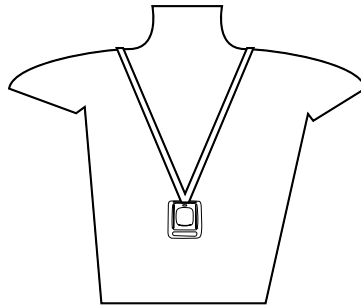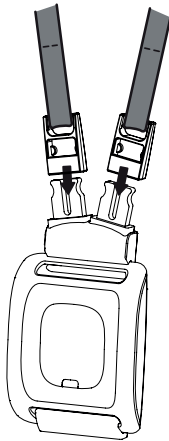

1

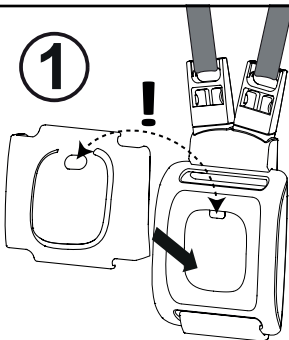


2

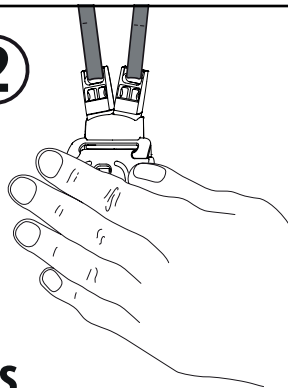


easy press

Designed to make the pendant easier to press, so you are only required to add on if you struggle to press the button.


①

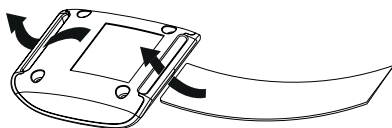


②

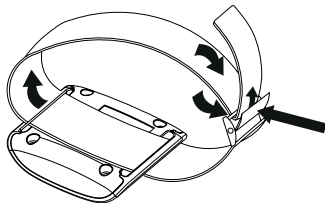

③

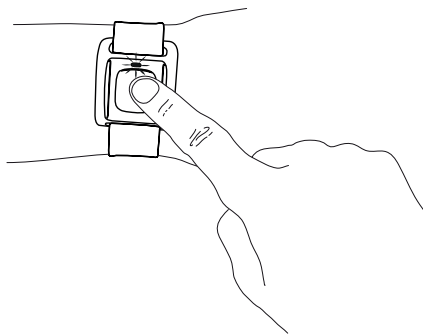


+

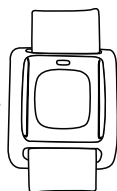

=

OK

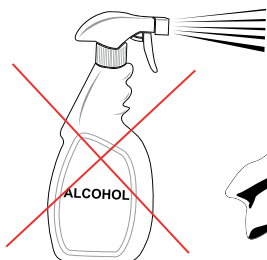

JUNE 2023

Careium
Aspinall House,
Walker Office Park,
Blackburn,
Lancashire,
BB1 2QE

Call: 0300 333 6511
Visit: www.careium.co.uk
Email: uk.info@careium.com

